

FILT'RAYScan

MANAGE AND SCAN POU FILT'RAY 2G FILTERS

Android & iOS

User Guide



FILT'RAYScan Traceability Management Tool

INTUITIVE MOBILE APPLICATION DEDICATED TO THE MANAGEMENT AND MAINTENANCE OF TERMINAL WATER FILTERS FILT'RAY 2G

- Manage and Maintain FILT'RAY 2G shower and terminal tap water filters in the field from any
 mobile device with an internet connection:
 - Android 5.0 or higher operating system
 (example: Samsung Xcover 4S is under Android 9).
 - iPhone and iPad version iOS 10 (example: iPhone 11 is on OS 13).
- Secure identification platform.





Using FILTRAYScan Features

- **Download & install FILTRAYSCAN** app on Android & iOS Devices
- Setup the device for automatic updates
- Login into FILT'RAYScan
- Dashboard Menu Overview
- Main Menu Overview
- Building Management
 Create Buildings & Floors / Create Bedrooms & Water Points
- Add Filters to Water points

 Scan filter barre code label and associate it with the Water Point
- Set Change out reminder alerts
- Create a pick list & organize the change out for the maintenance agents.
- Generate Reports

Save and send intervention export reports in PDF and Excel format.





Download & Install FILT'RAYScan on Android Devices



Click the **Play Store application** on your phone.


- Click the **Applications search bar**
- Enter "aqua-tools"
- Select "FILT'RAYScan"
- Click on "Install"

Follow the onscreen instructions to complete the installation.

Installation of the "FILT'RAYScan" application only takes for a few seconds.

Direct download link:

https://play.google.com/store/apps/details?id=com.aquatools.filtray&hl=en

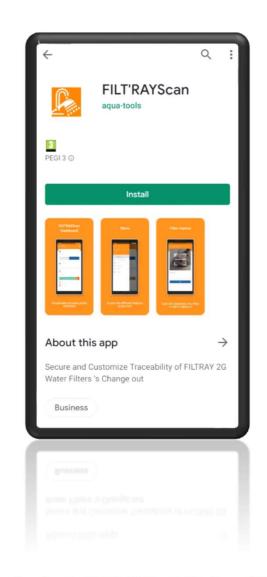





Download & Install FILT'RAYScan on iOS Devices



Open the **App Store** on your phone.


- Click the **App search bar**
- Enter "aqua-tools"
- Select "FILT'RAYScan"
- Click on "**Install**"

Follow the onscreen instructions to complete the installation.


Installation of the "FILT'RAYScan" application only takes for a few seconds.

Direct download link:

https://apps.apple.com/fr/app/filtrayscan/id1450244512?l=en

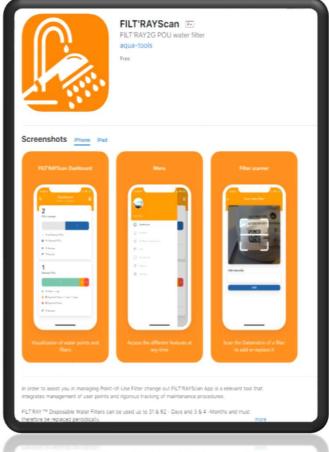



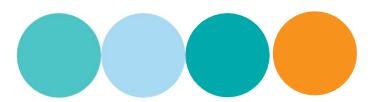

FILT'RAYScan automatic update on Devices

Android update:

Open the **Google Play Store** application Tap Settings **Menu**.

Touch Auto Update Applications.

Select the option:

"On any network to update applications using a Wi-Fi connection or mobile data".


Apple update:

On an iPhone, iPad, go to **Settings** then **[your name]**,


iTunes Store and App Store,

Then activate the Automatic updates option.



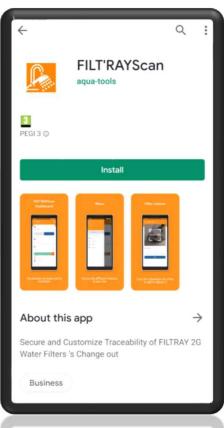




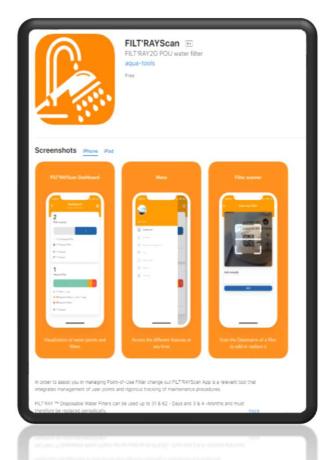




FILT'RAYScan Login Access

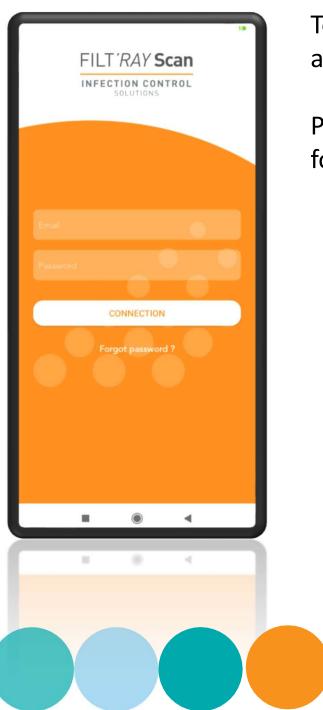

To get acces to the app enter your login information and select a profile.


Please contact your aqua-tools sales representative for login credentials.

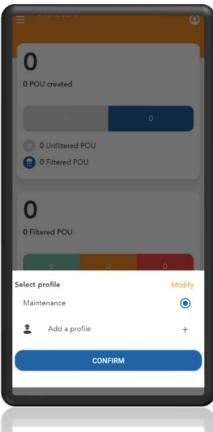





FILT'RAYScan Dashboard

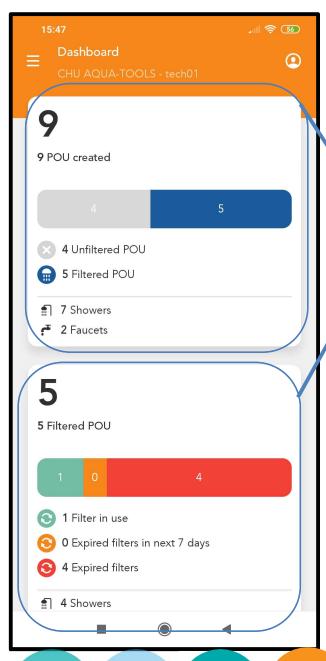

Dashboard Menu allows a global view of all the point-of-use water connections that have been created in the account.

Top window indicates how many points-of-use are created, what types of filters they require, how many have filters installed, and how many are running without filters.

Bottom window indicates duration time of installed point-of use water filters.

Note: The Dashboard is interactive, you can click on the bars to get the location of the points-of-use.

Colors codes:

Blue: Number of points-of-use with Water Filters

Grey: Water Points without water filters

Green: Number of points-of-use water filters RUNNING Orange: Number of filters NEEDING to be CHANGED soon


Red: Number of POUs with EXPIRED water filters



FILT'RAYScan Main Menu

To Open the Main Menu, click this icon on the top left:



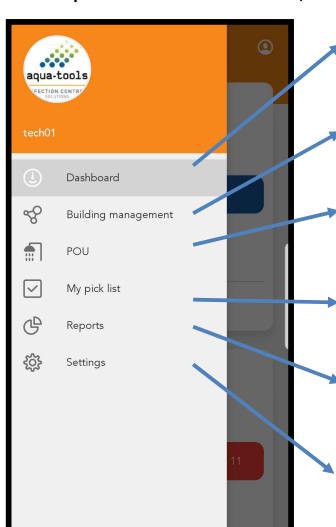

Dashboard gives a full overview of all the point-of-use connections registered on the account, the types of filters (shower/faucet), and their running times.


Building Management permits creation of all location layers (floors, departments, rooms), and matches the POU filters to the locations.

POU shows the FULL list of the points-of-use where water filters are running, provides a search/scope tool, allows activativation/deactivation of points-of-use, allows pick list creation before a water filter change-out tour.

My Pick List lists the filter types and locations for your the next filter change-out tour

Reports generates reports of change-outs to download in PDF or EXCEL format with results sorted by type of water filters, locations, and reasons.

Settings allows setup of change-out reminder alerts, notifications, account data updates, and sign out.

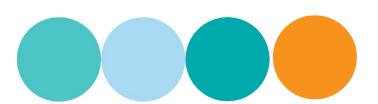






Building Management

Create Buildings and Floors


To Open the Main Menu, find on top this icon

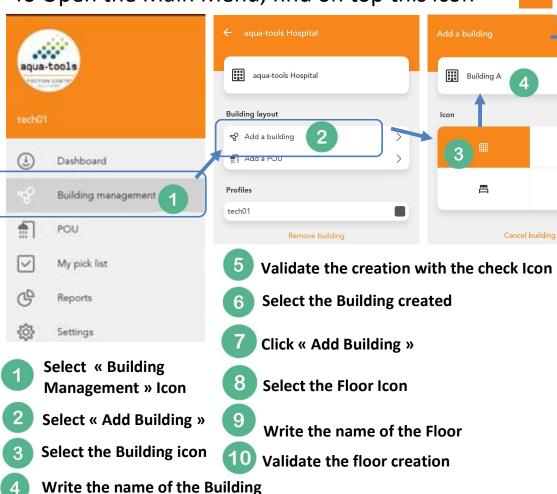

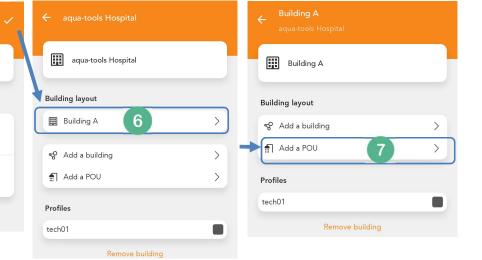

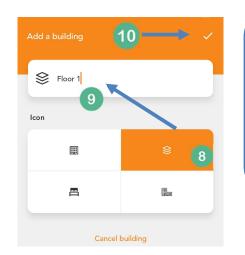



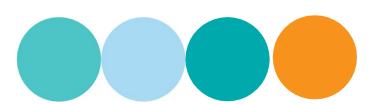




Building Management

Create Room Locations & Water Points in Each Room


To Open the Main Menu, find on top this icon




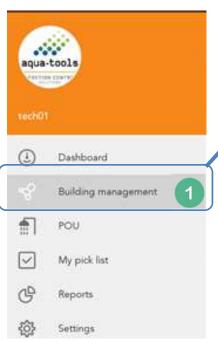





- 2 Select the Building
- 3 Select the Floor
- Slect Add a Building
- Select the Bedroom Icon
- 6 Write the name of the Bedroom
- Validate the creation with the check Icon

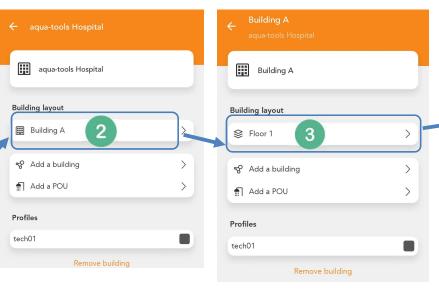

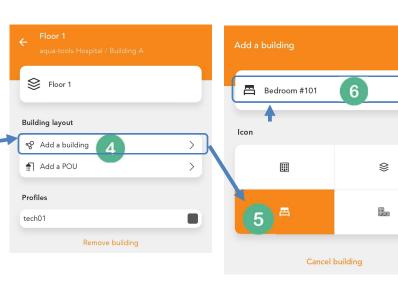

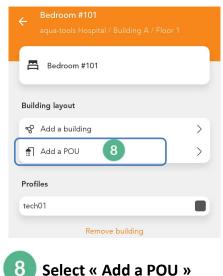

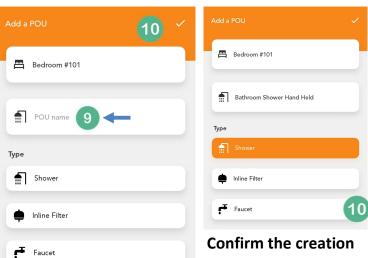








Add Filters to Water points

To Open the Main Menu, find on top this icon



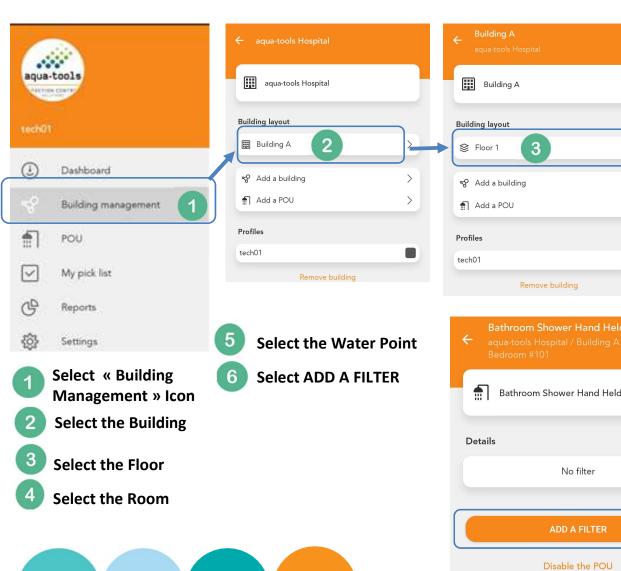

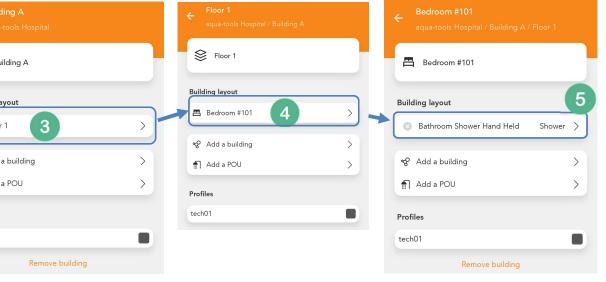


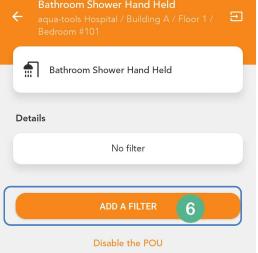





Building Management

Add Filters to Water points

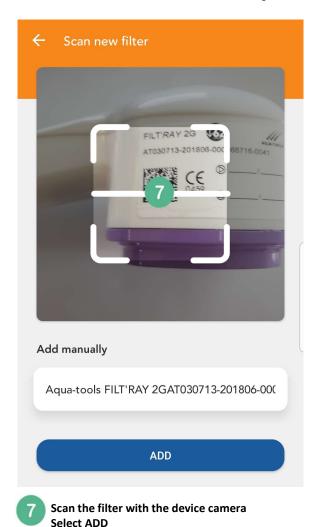

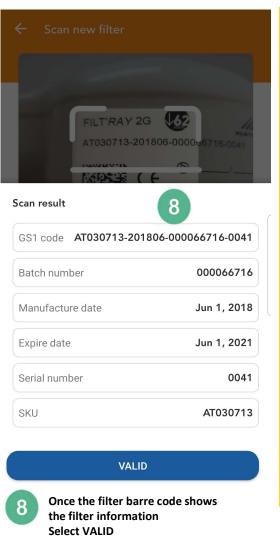

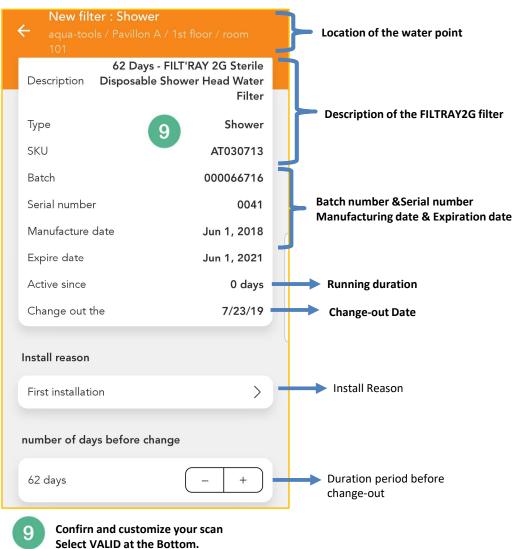


Note: The duration of the filters and the reason for the change-out can be customized at the step 9.







Setup Change-out reminder alerts


Open the Main Menu on top left of the app with this icon

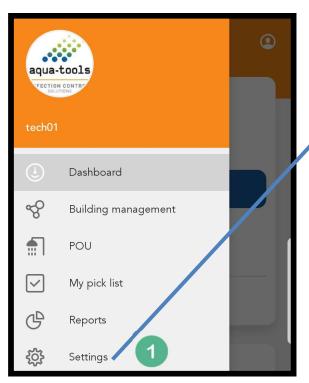

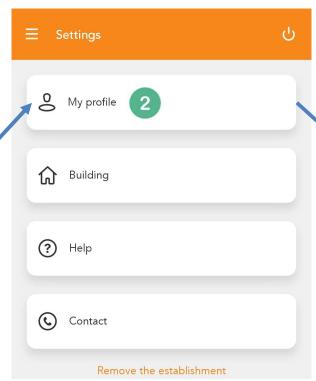

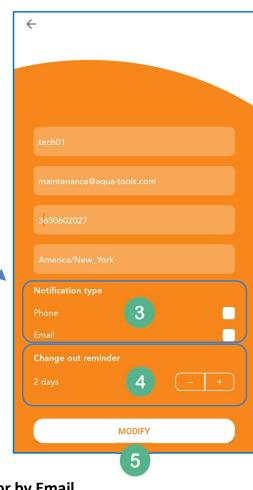

- Open the Main MENU and select Settings
- 2 Select « My PROFILE»

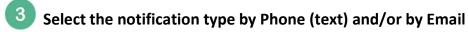

- 4 Set the number of days to be remided before the next change-out by + icons
- 5 Confirm by « MODIFY »






FILT'RAYScan POU Menu

View of all the filter point-of-use locations recorded in the account

- Find POUs by the search scope tool
- Activate/desactivate POU water filters
- Prepare the pick list before a change out

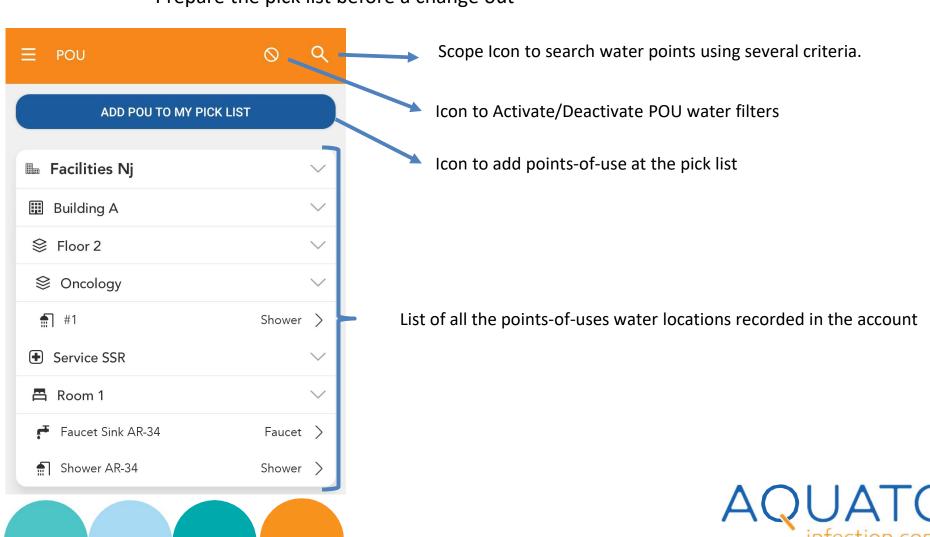






POU MENU: Create Pick List

Open the Main Menu on top of the app with this icon



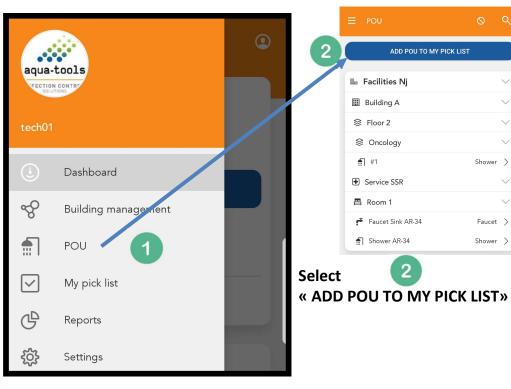

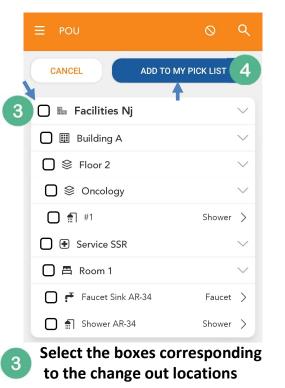

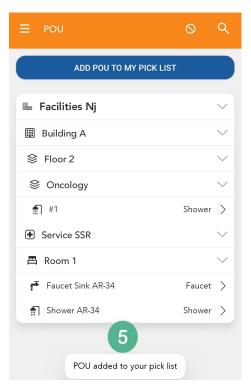

Open the Main MENU and select POU

- Confirm by selecting
 « ADD TO MY PICK LIST »
 - A Popup Message confirms POU added to your pick list








My Pick List: Change out

Open the Main Menu on top of the app with this icon




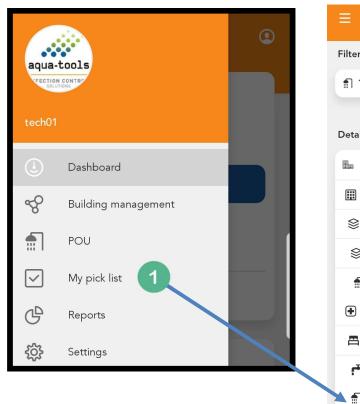

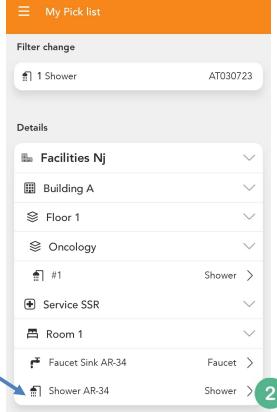


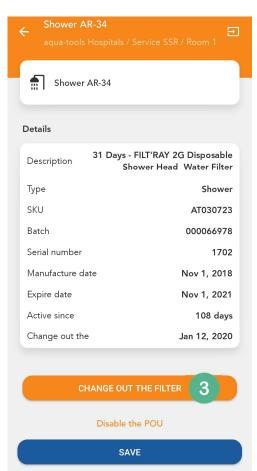



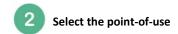




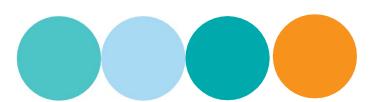




My Pick List: Change out



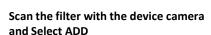

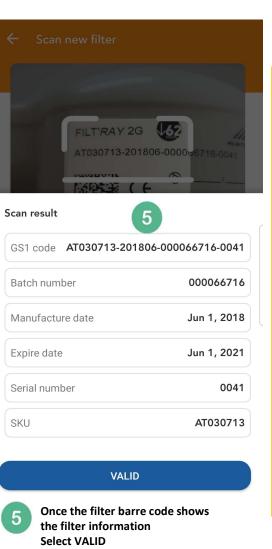

62 Days - FILT'RAY 2G Sterile Disposable Shower Head Water Description Filter Description of the FILTRAY2G filter Shower Type SKU AT030713 6 000066716 Batch Batch number & Serial number Serial number 0041 Manufacturing date & Expiration date Manufacture date Jun 1, 2018 Expire date Jun 1, 2021 Running duration Active since 0 days Change out the Change-out Date 7/23/19 Install reason Install Reason First installation number of days before change 62 days Duration to use product before Change-out Confirn and customize your scan

Select VALID at the Bottom.

Note: The duration of the product use and the reason for the change-out can be personalized at the step 9.







Filter Change-out Reports

Open the Main Menu on top of the app with this icon



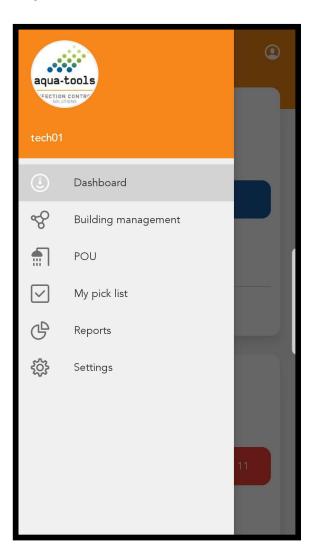

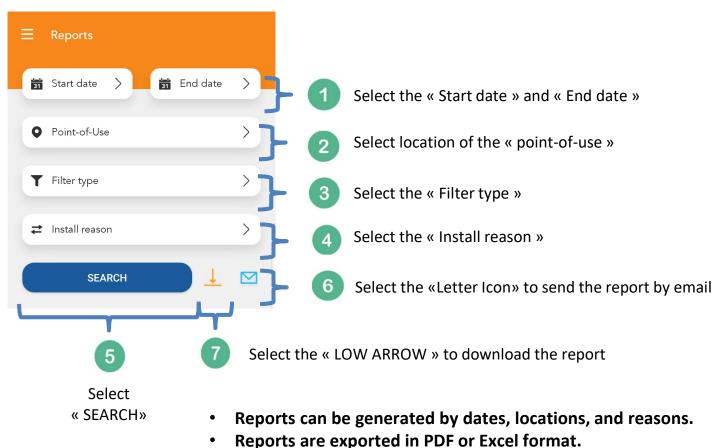

device





Reports can be sent to an email address and/or downloaded into the

FILT'RAYScan Help Menu

Open the Main Menu on top of the app with this icon





